

PLAMON

Optimal Plasma Device for Aesthetic Experts




Specialized for the Treatment of

- * Acne Care
- * Facial Lifting for Deep Wrinkles
- * Cosmetic Absorbent Functioning
- * Plasma Arc

Specifications

- Output Frequency : 250kHz
- Power Level : Step 1~9 (Analog wheel type)
- Output Mode : Continuous, Pulse
- Pulse Mode Frequency : 500Hz
- Power : Internal Li-Polymer Rechargeable Battery 3.7V, 500mAh
- This product is guaranteed for 1 year.
- Life Time : Full charged, Max 2 hours

Power Control by Analog Wheel

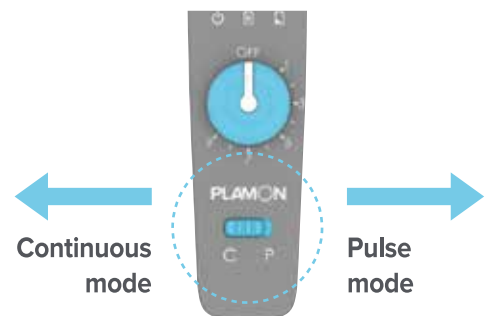


Support pulse mode of 500Hz as well as continuous mode.

Optimal Plasma Device for Aesthetic Experts

PLAMON

Treatments

Plasma Acne

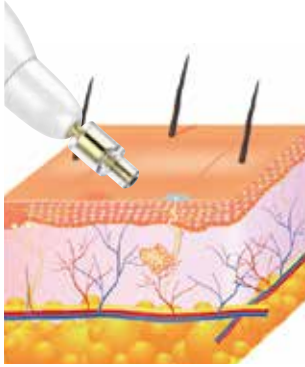


Acne Tip + Dedicated Film

When attach to the skin, a ring-shaped plasma forms at the tip of the film attached and releases adequate amount of radicals, ultraviolet (UV) and heat to help improve acne.

Acne Tip + Acne Patch


Plasma ring generates proper amount of radicals and UV rays under the acne patch.

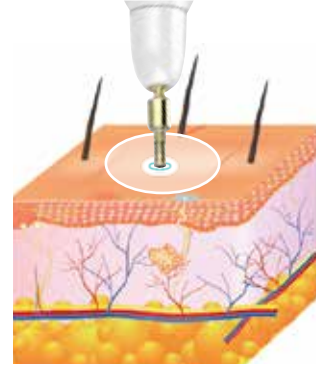

Plasma Lifting

Lifting Tip

- Support Pulse mode of 500Hz as well as Continuous mode.
- Plasma treatment with Pulse mode which is optimized for skin lifting collagen synthesis and wrinkles

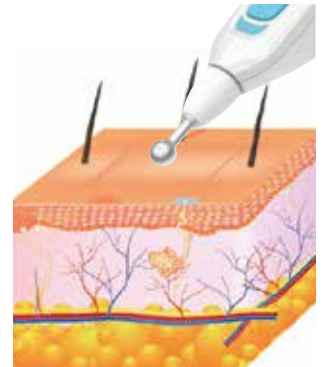

Cosmetic Absorbent

Fractional Tip & Permeating Tip

- 40kHz Continuous mode permeating tips with electroporation effect
- Fractional tips which are used together with Permeating tips maximize permeating effect

